


SIEMENS

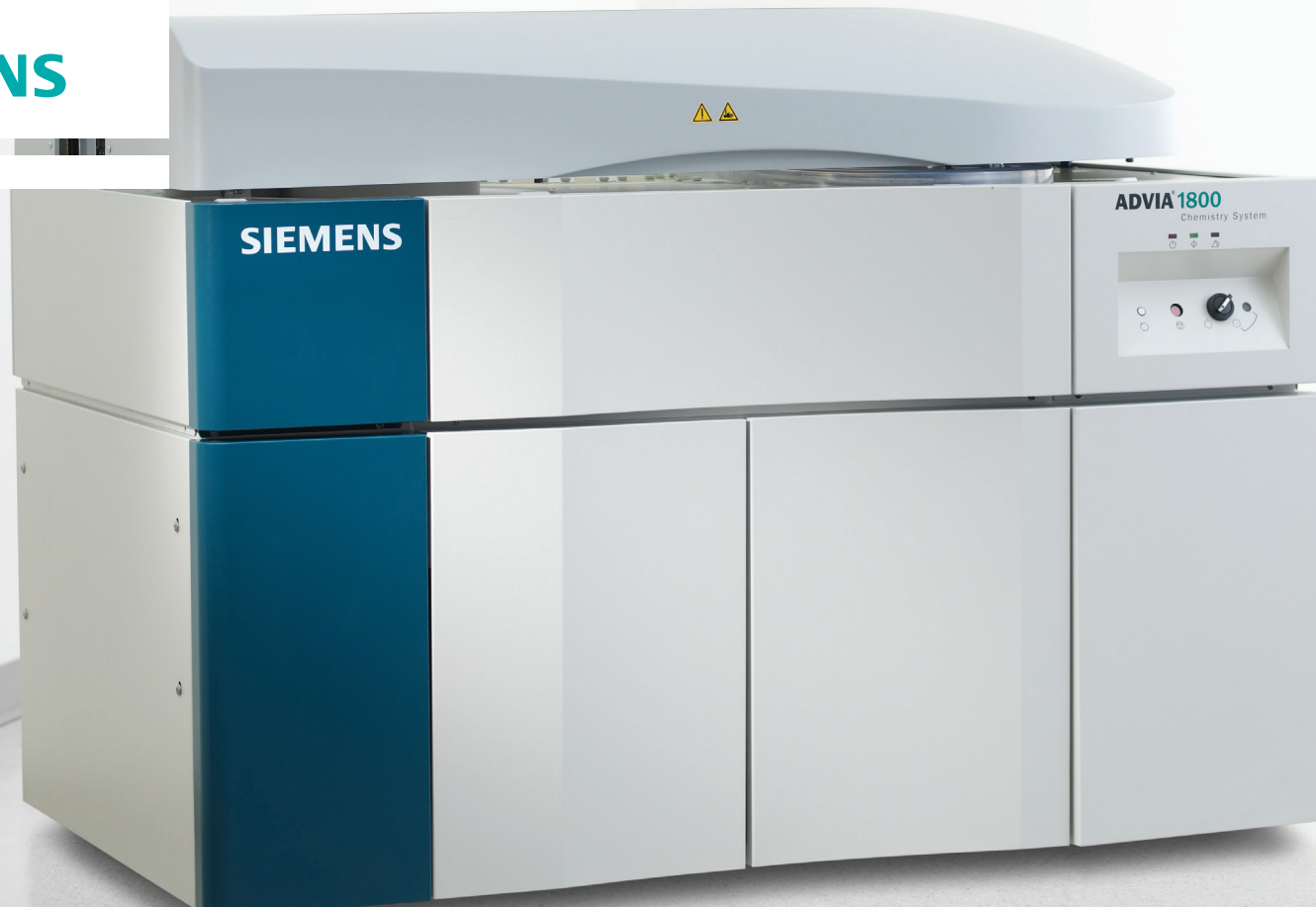

www.siemens.com/diagnostics

ADVIA Chemistry 1800 Assays Overview

Answers for life.

ADVIA® Chemistry 1800 Assays Overview

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
General Chemistry	Albumin BCG [ALB]	07622536 [10311832]	2220 [4 x 555]	S, P (Li Hep)	1–6 g/dL [1–6 g/dL]	BCG dye binding	60 days	Chemistry Cal 10312279	60 days	At calibration	10493936_EN Rev. D, 2011–01
General Chemistry	Albumin BCG [ALB_c]	04992820 [10335870]	6800 [4 x 1700]	S, P (Li Hep)	1–6 g/dL [1–6 g/dL]	BCG dye binding	60 days	Chemistry Cal 10312279	60 days	At calibration	10493937_EN Rev. C, 2011–01
General Chemistry	Albumin BCP [ALBP]	10711840 [10711840]	3200 [4 x 800]	S, P (EDTA, Li Hep)	0.6–8.0 g/dL [6–80 g/L]	BCG dye binding	60 days	Albumin BCP Cal 10711842	60 days	35 days	10699220_EN Rev. B, 2014-01
		10711841 [10711841]	6100 [4 x 1525]								
General Chemistry	Alkaline Phosphatase [ALPAMP]	03035814 [10309049]	2450 [7 x 350]	S, P (Li Hep)	0–1100 U/L [0–1100 U/L]	Modified IFCC	9 days with RCI [4 days w/o RCI]	Fixed factor	NA	Daily	10493939_EN Rev. F, 2014-06
General Chemistry	Alkaline Phosphatase [ALPA_c]	06860450 [10283340]	7440 [6 x 1240]	S, P (Li Hep)	5–1100 U/L [5–1100 U/L]	Modified IFCC	10 days with RCI [5 days w/o RCI]	Fixed factor	NA	Daily	10493938_EN Rev. B, 2011-02
General Chemistry	Alkaline Phosphatase [ALPDEA]	03036535 [10341131]	2485 [7 x 355]	S, P (Li Hep)	0–1600 U/L [0–1600 U/L]	DEA	40 days with RCI [8 days w/o RCI]	Fixed factor	NA	Daily	10493940_EN Rev. F, 2014-06 (OUS)
General Chemistry	Alanine Aminotransferase [ALT]	07501976 [10309500]	2520 [7 x 360]	S, P (Li Hep)	0–1100 U/L [0–1100 U/L]	Modified IFCC	60 days	Fixed factor	NA	Daily	10493941_EN Rev. F, 2014-04
		03036926 [10318168]	4020 [6 x 670]								
General Chemistry	Alanine Aminotransferase [ALTP5P]	07501976 [10309500]	2520 [7 x 360]	S, P (Li Hep)	5–1000 U/L [5–1000 U/L]	IFCC	30 days	Fixed factor	NA	Daily	10493944_EN Rev. G, 2014-06
		03036926 [10318168]	4020 [6 x 670]								
General Chemistry	Alanine Aminotransferase [ALT_c]	06860469 [10283341]	5100 [6 x 850]	S, P (Li Hep)	8–1100 U/L [8–1100 U/L]	Modified IFCC	45 days	Fixed factor	NA	Daily	10493942_EN Rev. D, 2011-11
General Chemistry	Alanine Aminotransferase [ALTP_c]	06860469 [10283341]	5100 [6 x 850]	S, P (Li Hep)	6–1000 U/L [6–1000 U/L]	Modified IFCC	30 days	Fixed factor	NA	Daily	10493943_EN Rev. D, 2011-01
General Chemistry	Ammonia [AMM]	00735335 [10286035]	200 [2 x 100]	S, P (EDTA, Li Hep)	10.0–1300 µg/dL [5.9–767.0 µmol/L]	Enzymatic	30 days	Tox/Ammonia Cal 10309217	30 days	Daily	10493945_EN Rev. G, 2014-04
General Chemistry	Amylase [AMYLAS]	07498401 [10309493]	1225 [7 x 175]	S, P (EDTA, Li Hep) U	0–1500 U/L [0–1500 U/L]	Ethylidene blocked-pNPG7	50 days	Fixed factor	NA	Daily	10493947_EN Rev. E, 2011-02
		03031177 [10341130]	2485 [7 x 355]								

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
General Chemistry	Aspartate Aminotransferase [AST]	07499718 [10309496]	2520 [7 x 360]	S, P (Li Hep)	0–1000 U/L [0–1000 U/L]	Modified IFCC	60 days	Fixed factor	NA	Daily	10493951_EN Rev. E, 2011-02
		03039631 [10341132]	4020 [6 x 670]								
General Chemistry	Aspartate Aminotransferase [ASTP5P]	03039631 10341132	2520 [7 x 360]	S, P (Li Hep)	5–1000 U/L [5–1000 U/L]	IFCC	30 days	Fixed factor	NA	Daily	10493954_EN Rev. F, 2012-04
		07499718 [10309496]	4020 [6 x 670]								
General Chemistry	Aspartate Aminotransferase [AST_c]	06860485 [10283343]	5100 [6 x 850]	S, P (Li Hep)	7–1000 U/L [7–1000 U/L]	Modified IFCC	60 days	Fixed factor	NA	Daily	10493952_EN Rev. D, 2011-01
General Chemistry	Aspartate Aminotransferase [ASTP_c]	06860485 [10283343]	5100 [6 x 850]	S, P (Li Hep)	8–1000 U/L [8–1000 U/L]	Modified IFCC	30 days	Fixed factor	NA	Daily	10493953_EN Rev. D, 2011-02
General Chemistry	Calcium [CA]	07506552 [10309504]	2506 [7 x 358]	S, P (Li Hep) U	S, P: 1.0–15.0 mg/dL [0.25–3.75 mmol/L] U: 1–30 mg/dL [0.25–7.5 mmol/L]	CPC	30 days with RCI [14 days w/o RCI]	Chemistry Cal 10312279	4 days with RCI (Daily w/o RCI)	At calibration	10493961_EN Rev. E, 2011-02
		04802312 [10286036]	4050 [6 x 675]								
General Chemistry	Calcium [CA_c]	06502821 [10282576]	5800 [4 x 1450]	S, P (Li Hep) U	S, P: 1.0–15.0 mg/dL [0.25–3.75 mmol/L] U: 1–30 mg/dL [0.25–7.5 mmol/L]	CPC	30 days with RCI [14 days w/o RCI]	Chemistry Cal 10312279	4 days with RCI (Daily w/o RCI)	At calibration	10493964_EN Rev. B, 2011-02
General Chemistry	Calcium [CA_2]	02189915 [10341116]	2030 [7 x 290]	S, P (Li Hep) U	S, P: 1.0–16.0 mg/dL [0.25–4.0 mmol/L] U: 1–32 mg/dL [0.25–8.0]	Arsenazo III	30 days	Chemistry Cal 10312279	24 days	At calibration	10493962_EN Rev. E, 2011-02
		02189699 [10316612]	4400 [8 x 550]								
General Chemistry	Calcium [CA_2c]	06860493 [10283344]	6800 [4 x 1700]	S, P (Li Hep) U	S, P: 1.0–16.0 mg/dL [0.25–4.0 mmol/L] U: 1–32 mg/dL [0.25–8.0]	Arsenazo III	42 days	Chemistry Cal 10312279	21 days	At calibration	10493963_EN Rev. B, 2011-02
General Chemistry	Carbon Dioxide [CO2_L]	07506781 [10326480]	1200 [4 x 300]	S, P (Li Hep)	10–40 mmol/L [1–40 mEq/L]	Enzymatic	30 days	CO2 Cal 10326423	Daily	Daily	10493973_EN Rev. F, 2014-04
		07987100 [10327377]	2400 [4 x 600]								
General Chemistry	Carbon Dioxide [CO2_c]	06860507 [10283345]	7520 [4 x 1880]	S, P (Li Hep)	1–40 mEq/L [1–40 mEq/L]	Enzymatic	30 days	CO2 Cal 10326423	Daily	Daily	10493972_EN Rev. A, 2010-10
General Chemistry	Cholesterol [CHOL_2]	10376501 [10376501]	2450 [7 x 350]	S, P (Li Hep)	10–675 mg/dL [0.26–17.48 mmol/L]	Enzymatic	60 days	Chemistry Cal 10312279	60 days	7 days	10493968_EN Rev. B, 2011-02
General Chemistry	Cholesterol [CHOL_c]	04993681 10335874]	8400 [4 x 2100]	S, P (Li Hep)	10–675 mg/dL [0.26–17.48 mmol/L]	Enzymatic	60 days	Chemistry Cal 10312279	60 days	7 days	10493969_EN Rev. B, 2011-02

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
General Chemistry	HDL Cholesterol [D-HDL]	07511947 [10309507]	1304 [4 x 326]	S, P (EDTA, Li Hep)	5–115 mg/dL [0.1–3.0 mmol/L]	Elimination/catalase	60 days	HDL/LDL Cal 10311918	30 days	At calibration	10493983_EN Rev. E, 2013-09
		06521559 [10283089]	4020 [6 x 670]								
General Chemistry	LDL Cholesterol [D LDL]	09796248 [10311891]	752 [4 x 188]	S, P (EDTA, Li Hep)	0–1000 mg/dL [0–25.9 mmol/L]	Elimination/catalase	14 days	HDL/LDL Cal 10311918	14 days	At calibration	10493981_EN Rev. D, 2011-02
General Chemistry	Cholinesterase [CHE]	01410545 [10315176]	480 [4 x 120]	S, P (Li Hep)	1,500–30,000 U/L [1,500–30,000 U/L]	Butrylthiocholine/5- thio-2-nitrobenzoate	7 days	Chemistry Cal 10312279	NA	Daily	10493966_EN Rev. D, 2011-02
General Chemistry	Creatine Kinase [CK NAC]	07498541 [10309494]	980 [7 x 140]	S, P (Li Hep)	0–1300 U/L [0–1300 U/L]	NAC activated, IFCC	10 days	Fixed factor	NA	Daily	10493970_EN Rev. D, 2011-02
		02096577 [10341111]	1080 [4 x 270]								
General Chemistry	Creatinine [CREA_2]	03039070 [10309050]	4020 [6 x 670]	S, P (Li Hep) U	S, P: 0.1–25.0 mg/dL [8.84–2210 µmol/L] U: 3.0–300 mg/dL [265–26,250 µmol/L]	Jaffe, alkaline picrate, kinetic with blank rate correction	20 days with RCI [9 days w/o RCI]	Chemistry Cal 10312279	3 days	At calibration	10493978_EN Rev. J, 2014-06
General Chemistry	Creatinine [CRE_2c]	06860515 [10283346]	9440 [4 x 2360]	S, P (Li Hep) U	S, P: 0.1–25.0 mg/dL [8.84–2210 µmol/L] U: 3.0–300 mg/dL [265–26,250 µmol/L]	Jaffe, alkaline picrate, kinetic with blank rate correction	35 days with RCI [15 days w/o RCI]	Chemistry Cal 10312279	20 days (3 days w/o RCI)	At calibration	10493976_EN Rev. D, 2013-01
General Chemistry	Enzymatic Creatinine [ECRE_2]	04992596 [10335869]	1400 [4 x 350]	S, P (EDTA, Li Hep) U	S, P: 0.1–30.0 mg/dL [9–2652 µmol/L] U: 1.00–300 mg/dL [88–26,520 µmol/L]	Creatininase	60 days	Chemistry Cal 10312279	60 days	At calibration	10493986_EN Rev. F, 2011-02
General Chemistry	Cystatin C [CYSC]	04851534 [10361930]	400 [4 x 100]	S, P (EDTA, Li Hep)	0.10–8.3 mg/L [0.10–8.3 mg/L]	Latex enhanced immunoturbidimetric	60 days	Cystatin C Cal 10376487	60 days	At calibration	10493980_EN Rev. E, 2011-02
General Chemistry	Direct Bilirubin [DBIL_2]	02189494 [10316610]	520 [4 x 130]	S, P (Li Hep)	0.1–15.0 mg/dL [2–256 µmol/L]	Vanadate oxidation	30 days	Chemistry Cal 10312279	30 days	At calibration	10493982_EN Rev. F, 2011-11
		02188390 [10341114]	2688 [4 x 672]								
General Chemistry	Gamma-Glutamyl Transferase [GGT]	07498649 [10309495]	980 [7 x 140]	S, P (Li Hep)	0–1200 U/L [0–1200 U/L]	Modified IFCC	30 days	Fixed factor	NA	Daily	10493990_EN Rev. E, 2014-05
		02011954 [10316298]	2700 [4 x 675]								
General Chemistry	Glucose [GLUH_3]	05001429 [10335891]	3960 [6 x 660]	S, P (EDTA, Li Hep) U,CSF	4–700 mg/dL [0.2–38.9 mmol/L]	Hexokinase	60 days	Chemistry Cal 10312279	60 days	60 days	10493992_EN Rev. D, 2012-02 (US)
				S, P (EDTA, Li Hep, Na FL/K oxalate) U,CSF							10698787_EN Rev. D, 2012-02 (OUS)

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
General Chemistry	Glucose [GLUH_c]	10377502 [10377502]	9400 [4 x 2350]	S, P (EDTA, Li Hep) U,CSF	4–700 mg/dL [0.2–38.9 mmol/L]	Hexokinase	60 days	Chemistry Cal 10312279	60 days	30 days	10493993_EN Rev. C, 2012-02 (US)
				S, P (EDTA, Li Hep, Na FL/K oxalate) U,CSF							10698786_EN Rev. C, 2012-02 (OUS)
General Chemistry	Glucose [GLUO]	10492319 [10492319]	2450 [7 x 350]	S, P (Li Hep) U,CSF	6–750 mg/dL [0.3–41.6 mmol/L]	Glucose Oxidase Trinder	60 days	Chemistry Cal 10312279	60 days	7days	10493994_EN Rev. E, 2012-02 (US)
				S, P (EDTA, Li Hep, Na FL/K oxalate) U,CSF							10698788_EN Rev. E, 2012-02 (OUS)
General Chemistry	Glucose [GLUO_c]	04993223 [10335872]	8480 [4 x 2120]	S, P (Li Hep) U,CSF	6–750 mg/dL [0.3–41.6 mmol/L]	Glucose Oxidase Trinder	60 days	Chemistry Cal 10312279	60 days	7 days	10493995_EN Rev. B, 2012-02 (US)
											10698789_EN Rev. B, 2012-02 (OUS)
General Chemistry	Inorganic Phosphorus [IP]	07506099 [10309503]	2506 [7 x 358]	S, P (Li Hep) U	S, P: 0–20 mg/dL [0–6.46 mmol/L] U: 4.0–200 mg/dL [1.29–64.60 mmol/L]	Phosphomolybdate/ UV	50 days	Chemistry Cal 10312279	50 days	7 days	10494002_EN Rev. F, 2014-05
		03056927 [10318211]	4020 [6 x 670]								
General Chemistry	Inorganic Phosphorus [IP_c]	06860523 [10283347]	10,200 [6 x 1700]	S, P (Li Hep) U	S, P: 0–20 mg/dL [0–6.46 mmol/L] U: 4.0–200 mg/dL [1.29–64.60 mmol/L]	Phosphomolybdate/ UV	60 days	Chemistry Cal 10312279	60 days	7 days	10494003_EN Rev. D, 2014-05
General Chemistry	Iron [IRON_2]	02194838 [1034118]	1015 [7 x 145]	S, P (Li Hep)	2–1000 µg/dL [0.3–179 µmol/L]	Ferrozine	60 days	Chemistry Cal 10312279	60 days	At calibration	10494004_EN Rev. E, 2012-05
		10377510 [10377510]	2100 [6 x 350]								
General Chemistry	Lactate Dehydrogenase [LDPL]	07502999 [10309502]	1790 [5 x 358]	S, P (Li Hep)	12–1000 U/L [12–1000]	Pyruvate/NADH	50 days	Fixed factor	NA	Daily	10494008_EN Rev. E, 2012-04
		03030863 [10341129]	4020 [6 x 670]								
General Chemistry	Lactate Dehydrogenase [LDLP]	07502115 [10309501]	1918 [7 x 274]	S, P (Li Hep)	20–700 U/L [20–700 U/L]	Lactate/NAD	30 days	Fixed factor	NA	Daily	10494007_EN Rev. D, 2011-02
		03029628 [10341128]	4020 [6 x 670]								
General Chemistry	Lactate [LAC]	07109944 [10325776]	300 [3 x 100]	P (EDTA, Hep Citrate, Ox/Fl)	0.5–110.0 mg/dL [0.06–12.21 mmol/L]	Lactate oxidase, colorimetric	45 days	Special CHEM Cal 10311926	45 days	At calibration	10494006_EN Rev. F, 2014-06

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
General Chemistry	Lipase [LIP]	01984894 [10311896]	640 [4 x 160]	S, P (Li Hep)	0–700 U/L [0–700 U/L]	Colorimetric rate	10 days	Special CHEM Cal 10311926	10 days	At calibration	10494009_EN Rev. E, 2011-08
General Chemistry	Magnesium [MG]	07508865 [10309506]	980 [7 x 140]	S, P (Li Hep) U	S, P: 0.7–5.0 mg/dL [0.29–2.06 mmol/L] U: 0.9–25.0 mg/dL [0.37–10.28]	Xyldyl blue	7 days with RCI [4 days w/o RCI]	Chemistry Cal 10312279	4 days	At calibration	10494014_EN Rev. G, 2014-06
		02271328 [10341119]	1120 [4 x 280]								
General Chemistry	Magnesium [MG_c]	04992960 [10335871]	7620 [6 x 1270]	S, P (Li Hep) U	S, P: 0.7–5.0 mg/dL [0.29–2.06 mmol/L] U: 0.9–25.0 mg/dL [0.37–10.28]	Xyldyl blue	10 days with RCI [4 days w/o RCI]	Chemistry Cal 10312279	10 days with RCI (4 days w/o RCI)	At calibration	10494015_EN Rev. D, 2014-06
General Chemistry	Pancreatic Amylase [PAMY]	01410820 [10315178]	910 [7 x 130]	S	0–1500 U/L [0–1500]	Antibody-blocked colorimetric amylase	60 days	Fixed factor	NA	Daily	10494015_EN Rev. D, 2014-06
General Chemistry	Total Iron Binding Capacity [TIBC]	03940010 [10309071]	400 [4 x 100]	S	40–670 µg/dL [7.2–119.9 umol/L]	Sequential release and uptake of iron	7 days with RCI	Special CHEM Cal 10311926	7 days with RCI	At calibration	10494041_EN Rev. D, 2011-02
General Chemistry	Total Bilirubin [TBIL_2]	02189346 [10341115]	1424 [4 x 356]	S, P (Li Hep)	0.1–35.0 mg/dL [2–598 µmol/L]	Vanadate oxidation	60 days	Chemistry Cal 10312279	60 days	At calibration	10494037_EN Rev. E, 2011-02
		02188064 [10341113]	2688 [4 x 672]								
General Chemistry	Total Protein [TP]	01803350 [10311878]	3400 [4 x 850]	S, P (Li Hep)	2.0–12.0 g/dL [20–120 g/L]	Biuret	60 days with RCI [21 days w/o RCI]	Chemistry Cal 10312279	30 days with RCI (21 days w/o RCI)	At calibration	10494044_EN Rev. H, 2014-06
General Chemistry	Total Protein [TP_c]	06860531 [10283348]	12,240 [4 x 3060]	S, P (Li Hep)	2.0–12.0 g/dL [20–120 g/L]	Biuret	60 days with RCI [30 days w/o RCI]	Chemistry Cal 10312279	30 days with RCI (21 days w/o RCI)	At calibration	10494045_EN Rev. C, 2011-01
General Chemistry	Triglycerides [TRIG_2]	10335892 [10335892]	1432 [4 x 358]	S, P (EDTA)	10–550 mg/dL [0.11–6.22 mmol/L]	GPO, Trinder without serum blank	60 days	Chemistry Cal 10312279	60 days	14 days	10494049_EN Rev. B, 2014-07
General Chemistry	Triglycerides [TRIG_c]	10697575 [10697575]	4252 [4 x 1063]	S, P (EDTA)	10–550 mg/dL [0.11–6.22 mmol/L]	GPO, Trinder without serum blank	60 days	Chemistry Cal 10312279	60 days	14 days	10494050_EN Rev. B, 2014-07
General Chemistry	Urea Nitrogen [UN]	03040257 [10309051]	4020 [6 x 670]	S, P (Li Hep)	5–150 mg/dL [1.8–53.6 mmol/L]	Urease with GLDH	30 days	Chemistry Cal 10312279	30 days	At calibration	10494054_EN Rev. D, 2011-01
General Chemistry	Urea Nitrogen [UN_c]	06860558 [10283349]	9440 [4 x 2360]	S, P (Li Hep)	5–150 mg/dL [1.8–53.6 mmol/L]	Urease with GLDH	30 days	Chemistry Cal 10312279	30 days	7 days	10494055_EN Rev. B, 2011-02
General Chemistry	Uric Acid [UA]	07497014 [10309492]	980 [7 x 140]	S, P (Li Hep)	0.5–20.0 mg/dL [30–1190 umol/L]	Uricase/Peroxidase	60 days	Chemistry Cal 10312279	60 days	7 days	10494051_EN Rev. E, 2011-02
		03051305 [10341133]	4020 [6 x 670]								
General Chemistry	Uric Acid [UA_c]	04993452 [10335873]	7500 [6 x 1250]	S, P (Li Hep)	0.5–20.0 mg/dL [30–1190 umol/L]	Uricase/Peroxidase	60 days	Chemistry Cal 10312279	60 days	7 days	10494052_EN Rev. C, 2014-06
General Chemistry	Total Protein [UPRO_2]	05000171 [10335885]	776 4 x 194	U CSF	U: 1.1–125 mg/dL [11–1250 mg/L] CSF: 1.2–125 mg/dL [12–1250 mg/L]	Dye binding	60 days with RCI [49 days w/o RCI]	Urine Protein Cal 10311930	30 days with RCI (14 days w/o RCI)	At calibration	10699690_EN Rev. E, 2014-06

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
Diabetes	Fructosamine [FRUCT]	04862501 [10361941]	200 [2 x 100]	S, P (EDTA, Li Hep)	30–1000 µmol/L [30–1000 µmol/L]	Enzymatic reaction	60 days	Fructosamine Cal 10285684	28 days	At calibration	10697377_EN Rev. B, 2015-01
Diabetes	Hemoglobin A1c [A1c_3]	10379673 [10379673]	160 [2 x 80]	S, P (EDTA, Li Hep)	A1c 1.0–8.83 µmol/L tHb 7–24 g/dL HbA1c 8–144 mmol/mol	Competition between hydrolysed sample HbA1c and agglutinator for anti_HbA1c antibody	30 days with RCI [7 days w/o RCI]	A1c_3 Cal 10491408	14 days with RCI (7 days w/o RCI)	At calibration	10493997_EN Rev. G, 2013-09
		10485591 [10485591]	800 [4 x 200]								
Diabetes	Hemoglobin A1c [A1c_3M]	10379673 [10379673]	160 [2 x 80]	S, P (EDTA, Li Hep)	A1c 1.0–8.83 µmol/L tHb 7–24 g/dL HbA1c 8–144 mmol/mol	Competition between hydrolysed sample HbA1c and agglutinator for anti_HbA1c antibody	30 days with RCI (7 days w/o RCI)	A1c_3 Cal 10491408	14 days with RCI (7 days w/o RCI)	At calibration	10485703_EN Rev. B, 2011-10
		10485591 [10485591]	800 [4 x 200]								
Diabetes	Microalbumin [uALB_2]	03051194 [10318197]	420 [4 x 105]	U	0.3–42.0 mg/dL [3–420] mg/L]	PEG–enhanced immunoturbidimetric	60 days	Microalbumin Cal 10337470	60 days	At calibration	10699695_EN Rev. F, 2014-06
Drugs of Abuse	Amphetamines [AMPH_2]	10378758 [10378758]	760 [4 x 190]	U	Semi-Quantitative Cutoff 300 ng/mL [100–1000 ng/mL] 500 ng/mL [150–1800 ng/mL] 1000 ng/mL [265–1800 ng/mL]	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit® Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10493946_EN Rev. C, 2011-11
Drugs of Abuse	Barbiturates [BARB_2]	10378759 [10378759]	760 [4 x 190]	U	Semi-Quantitative Cutoff 200 ng/mL [40–750 ng/mL] 300 ng/mL [40–750 ng/mL]	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10493957_EN Rev. C, 2012-01
Drugs of Abuse	Benzodiazepine [BENZ_2]	10378760 [10378760]	760 [4 x 190]	U	Semi-Quantitative Cutoff 200 ng/mL [50–900 ng/mL] 300 ng/mL [50–900 ng/mL] Qualitative	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10493958_EN Rev. C, 2011-11
Drugs of Abuse	Cannabinoid [THC_2]	10378766 [10378766]	760 [4 x 190]	U	Semi-Quantitative Cutoff 20 ng/mL [11–55 ng/mL] 50 ng/mL [15–113 ng/mL] 100 ng/mL [20–180 ng/mL] Qualitative	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10494039_EN Rev. C, 2011-11

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
Drugs of Abuse	Cocaine [COCA_2]	10378761 [10378761]	760 [4 x 190]	U	Semi-Quantitative Cutoff 150 ng/mL [36–900 ng/mL] 300 ng/mL [36–900 ng/mL]	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10493974_EN Rev. C, 2011-11
Drugs of Abuse	Ecstasy [XTC]	10378735 [10378735]	760 [4 x 190]	U	Semi-Quantitative Cutoff 300 ng/mL [75–500 ng/mL] 500 ng/mL [75–1000 ng/mL]	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit II Plus Ecstasy Cals/Controls Level 0 9A509UL Level 1 9X529UL Level 2 9X549UL Level 3 9X569UL Level 4 9X589UL	15 days	Daily	10494061_EN Rev. E, 2014-08
Drugs of Abuse	Methadone [MDN_2]	10378762 [10378762]	760 [4 x 190]	U	Semi-Quantitative Cutoff 150 ng/mL [40–450 ng/mL] 300 ng/mL [90–900 ng/mL] Qualitative	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10494011_EN Rev. C, 2011-11
Drugs of Abuse	Methadone Metabolite [METMTB]	05038055 [10335933]	1000 [4 x 250]	U	Cutoff 1000 ng/mL 0–1000 ng/mL [0–1000 ng/mL]	Homogenous competitive enzyme immunoassay	45 days	Methadone Metabolite Urine Cal 10283006 Negative Drug Urine Cal 10313977	14 days	Daily	10494013_EN Rev. D, 2011-02
Drugs of Abuse	Opiates [OP_2]	10378763 [10378763]	760 [4 x 190]	U	Semi-Quantitative Cutoff 300 ng/mL [75–1800 ng/mL] 2000 ng/mL [320–3800 ng/mL] Qualitative	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10494017_EN Rev. D, 2014-09
Drugs of Abuse	Phencyclidine [PCP_2]	10378764 [10378764]	760 [4 x 190]	U	Semi-Quantitative Cutoff 25 ng/mL [4–90 ng/mL] Qualitative	Enzyme multiplied immunoassay technique (EMIT)	30 DAYS	Emit Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10494017_EN Rev. D, 2014-09
Drugs of Abuse	Propoxyphene [PPX_2]	10378765 [10378765]	760 [4 x 190]	U	Semi-Quantitative Cutoff 300 ng/mL [103–450 ng/mL] Qualitative	Enzyme multiplied immunoassay technique (EMIT)	30 days	Emit Cals/Controls Level 0 9A509UL Level 1 9A529UL Level 2 9A549UL Level 3 9A569UL Level 4 9A589UL Level 5 9A609UL	30 days	Daily	10494025_EN Rev. C, 2011-11

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
ISE	Chloride [CL]	07097504 [10309439]		S, P, U	S, P: 15–200 mEq/L [15–200 mmol/L] U: 15–400 mEq/L [15–400 mmol/L]	Indirect potentiometric		ISE STD Serum 10311903 ISE STD Urine 10311904	Daily	NA	10493971_EN Rev. E, 2011-01
ISE	Potassium [K]	06135445 [10309440]		S, P, U	S, P: 1–10 mEq/L [1–10 mmol/L] U: 3–300 mEq/L [3–300 mmol/L]	Indirect potentiometric		ISE STD Serum 10311903 ISE STD Urine 10311904	Daily	NA	10494005_EN Rev. E, 2011-01
ISE	Sodium [NA]	03092699 [10309441]		S, P, U	S, P: 100–200 mEq/L [100–200 mmol/L] U: 10–400 mEq/L [10–400 mmol/L]	Indirect potentiometric		ISE STD Serum 10311903 ISE STD Urine 10311904	Daily	NA	10494016_EN Rev. E, 2011-01
Specific Proteins	Alpha-1 - Acid Glycoprotein [AAG]	04862129 [10361936]	200 [2 x 100]	S, P (EDTA, Li Hep)	5.0–675) mg/dL [0.05–6.75) g/L]	PEG immunturbidimetric	60 days	Liquid Specific Protein Cal 10311929	28 days with RCI (14 days w/o RCI)	At calibration	10493932_EN Rev. B, 2011-11
Specific Proteins	Alpha-1- Antitrypsin [AAT]	03059241 [10318221]	200 [2 x 100]	S, P (Li Hep)	5.0–540 mg/dL [0.05–5.40 g/L]	PEG immunturbidimetric	42 days	Liquid Specific Protein Cal 10311929	7 days	At calibration	10493933_EN Rev. E, 2011-02
Specific Proteins	Anti-Streptolysin O [ASO_2]	03055483 [10318208]	200 [2 x 100]	S, P (Li Hep)	25–1000 IU/mL [25–1000 IU/mL]	Latex-enhanced immunoturbidimetric	60 days	Liquid Specific Protein Cal 10311929	60 days	At calibration	10493950_EN Rev. D, 2011-02
Specific Proteins	Apolipoprotein A-1 [APO A1]	03050910 [10318196]	150 [2 x 75]	S, P (Li Hep)	15–260 mg/dL [0.15–2.6) g/L]	PEG immunturbidimetric	14 days	APO A1/B Cal 10311898	14 days	Daily	10493948_EN Rev. F, 2011-01
Specific Proteins	Apolipoprotein B [APO B]	03055068 [10318207]	170 [2 x 85]	S, P (Li Hep)	15–260) mg/dL [0.15–2.6) g/L]	PEG immunturbidimetric	30 days	APO A1/B Cal 10311898	30 days	Daily	10493949_EN Rev. F, 2011-01
Specific Proteins	Beta-2 Microglobulin [B2M]	04862161 [10361937]	200 [2 x 100]	S, P (EDTA, Li Hep)	0.25–18.00 mg/L [0.25–18.00 mg/L]	Latex-enhanced immunoturbidimetric	60 days	B2M Cal 10285680	21 days	Daily	10493955_EN Rev. B, 2012-05
Specific Proteins	C- Reactive Protein [CRP_2]	06522059 [10283129]	1000 [4 x 250]	S, P (Li Hep)	4–336 mg/L [0.4–33.6 mg/dL]	Latex-enhanced immunoturbidimetric	60 days	CRP_2 Cal 10337471	60 days	At calibration	10493979_EN Rev. E, 2014-06
Specific Proteins	Complement C3 [C3]	03060444 [10309053]	200 [2 x 100]	S	0.5–550 mg/dL [0.01–5.50 g/L]	PEG-enhanced immunoturbidimetric	42 days	Liquid Specific Protein Cal 10311929	14 days	At calibration	10493959_EN Rev. E, 2011-02
Specific Proteins	Complement C4 [C4]	03060770 [10309054]	200 [2 x 100]	S	0.4–150 mg/dL [0.01–1.50 g/L]	PEG-enhanced immunoturbidimetric	42 days	Liquid Specific Protein Cal 10311929	14 days	At calibration	10493960_EN Rev. D, 2011-02
Specific Proteins	Ferritin [FRT]	04862277 [10361938]	800 [4 x 200]	S, P (EDTA, Li Hep)	6.0–500.0 ng/mL [6.0–500) µg/L]	Latex microparticle- enhanced immunoturbidimetric	60 days	Liquid Specific Protein Cal 10311929	60 days	At calibration	10493988_EN Rev. C, 2013-08
		10494078 [10494078]	3200 [4 x 800]								
Specific Proteins	Haptoglobin [HAPT]	06149314 [10323850]	200 [4 x 50]	S	1–410 mg/dL [0.1–4.10 g/L]	PEG-enhanced immunoturbidimetric	30 days (180 test kit)	Liquid Specific Protein Cal 10311929	30 days (180 test kit)	At calibration	10493996_EN Rev. E, 2011-02
		0831008 [10314164]	540 [3 x 180]				7 days (50 test kit)		7 days (50 test kit)		

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
Specific Proteins	High Sensitivity C-Reactive Protein [hsCRP]	06837459 [10283282]	440 [2 x 220]	S, P (EDTA, Li Hep)	0.16–10.0 mg/L [0.16–10.0 mg/L]	Latex-enhanced immunoturbidimetric	60 days	hsCRP Cal 10335897	28 days	At calibration	10493998_EN Rev. D, 2014-05
Specific Proteins	Immunoglobulin A [IGA_2]	02194102 [10309033]	375 [5 x 75]	S, P (Li Hep)	33.0–660) mg/dL [0.33–6.60) g/L]	PEG-enhanced immunoturbidimetric	60 days	Liquid Specific Protein Cal 10311929	60 days	At calibration	10493999_EN Rev. D, 2011-02
Specific Proteins	Immunoglobulin G [IGG_2]	02193432 [10309031]	360 [4 x 90]	S, P (Li Hep)	140–4100 mg/dL [1.40–41.00 g/L]	PEG-enhanced immunoturbidimetric	60 days	Liquid Specific Protein Cal 10311929	60 days	At calibration	10494000_EN Rev. D, 2011-02
Specific Proteins	Immunoglobulin M [IGM_2]	02193483 10309032]	360 [4 x 90]	S, P (EDTA, Li Hep)	21.0–410 mg/dL [0.21–4.10 g/L]	PEG-enhanced immunoturbidimetric	60 days	Liquid Specific Protein Cal 10311929	60 days	At calibration	10494001_EN Rev. F, 2011-03
Specific Proteins	Lipoprotein a [LPA]	04862528 [10361942]	200 [2 x 100]	S, P (Li Hep)	10.00–85.00 mg/dL [0.100–0.850 g/L]	Latex microparticle- enhanced immunoturbidimetry	30 days	LPA Cal 10285017	30 days	At calibration	10697439_EN Rev. C, 2013-06
Specific Proteins	Myoglobin [MYO]	04862455 [10361939]	200 [2 x 100]	S, P (Li Hep)	22.0–680.0 ng/mL [22.0–680.0 ug/L]	Latex microparticle- enhanced immunoturbidimetry	60 days	Myoglobin Cal 10285682	21 days	At calibration	10697501_EN Rev. A, 2012-10
Specific Proteins	Prealbumin [PREALB]	03057680 [10318214]	200 [2 x 100]	S	5.0–80.0 mg/dL [50–800 mg/L]	PEG enhanced immunoturbidimetric	35 days	Liquid Specific Protein Cal 10311929	7 days	At calibration	10494026_EN Rev. E, 2011-02
Specific Proteins	Rheumatoid Factor [RF]	00655137 [10311894]	200 [2 x 100]	S	3.5–110 IU/mL [3.5–10 IU/mL]	Latex-enhanced immunoturbidimetric	30 days	Liquid Specific Protein Cal 10311929	14 days	At calibration	10494028_EN Rev. E, 2011-02
Specific Proteins	Transferrin [TRF]	03059160 [10318218]	440 [4 x 110]	S, P (EDTA, Li Hep)	0.01–550 mg/dL [0.0001–5.50) g/L]	PEG-enhanced immunoturbidimetric	30 days	Liquid Specific Protein Cal 10311929	14 days	At calibration	10494046_EN Rev. F, 2014-06
Specific Proteins	Wide Range C-Reactive Protein [wrCRP]	03108390 [10318309]	440 [2 x 220]	S, P (Li Hep)	0.012–16.4 mg/dL [0.12–164 mg/L]	Latex-enhanced immunoturbidimetric	35 days	WR-CRP Cal 10313256	21 days	At calibration	10494060_EN Rev. F, 2014-06
		0829585 [10341101]	2205 [7 x 315]								
TDM	Carbamazepine [CARB_2]	10376342 [10376342]	400 [4 x 100]	S, P (EDTA, Hep Citrate, Ox/FI)	1.0–22.2 µg/mL [4.2–93.9 µmol/L]	Enzyme multiplied immunoassay technique (EMIT)	30 days with RCI [5 days w/o RCI]	Drug Cal 2 10376771	30 days	Daily	10493965_EN Rev. B, 2011-02
TDM	Digoxin [DIG]	08085984 [10327549]	800 [4 x 200]	S, P (Li Hep)	0.1–5.0 ng/mL [0.13–6.40 nmol/L]	Latex-enhanced immunoturbidimetric	30 days	TDM Cal 10311922	7 days	At calibration	10493984_EN Rev. E, 2011-02
TDM	Gentamicin [GENT_2]	10377507 [10377507]	400 [4 x 100]	S, P (EDTA, Hep Citrate, Ox/FI)	0.4–13.6 µg/mL [0.9–29.4 µmol/L]	Enzyme multiplied immunoassay technique (EMIT)	30 days with RCI [14 days w/o RCI]	Drug Cal 2 10376771	30 days with RCI (14 days w/o RCI)	Daily	10493989_EN Rev. B, 2011-02
TDM	Lithium [LITH]	00141230 [10312850]	200 2 x 100	S, P (EDTA, Li Hep)	0.1–3.0 mmol/L [0.1–3.0 mmol/L]	Colorimetric	21 days with RCI [30 days w/o RCI]	Chemistry Cal 10312287	4 days	At calibration	10494010_EN Rev. F, 2014-06
TDM	Phenobarbital [PHNB_2]	10377505 [10377505]	400 [4 x 100]	S, P (EDTA, Hep Citrate, Ox/FI)	1.0–88.0) µg/mL [4.3–379.3 µmol/L]	Enzyme multiplied immunoassay technique (EMIT)	28 days	Drug Cal 1 10376770	28 days	Daily	10494023_EN Rev. B, 2011-02

Discipline [Disease State]	Assay Name [Test Mnemonic]	REF [SMN]	Tests [Kit Config.]	Specimen Types	Analytical Measuring Range [SI Units]	Assay Principle	On Board Stability	Calibrator	Calibration Interval	Reagent Blank (RBL) Frequency	IFU Reference
TDM	Phenytoin [PHNY_2]	10377506 [10377506]	400 [4 x 100]	S, P (EDTA, Hep Citrate, Ox/Fl)	1.0–40.0 µg/mL [4.0–158.4 µmol/L]	Enzyme multiplied immunoassay technique (EMIT)	28 days	Drug Cal 1 10376770	28 days	Daily	10494024_EN Rev. B, 2011-02
TDM	Theophylline [THEO_2]	10377503 [10377503]	400 [4 x 100]	S, P (EDTA, Hep Citrate, Ox/Fl)	1.5–40.0 µg/mL [8.3–222.0 µmol/L]	Enzyme multiplied immunoassay technique (EMIT)	28 days with RCI [7 days w/o RCI]	Drug Cal 1 10376770	28 day with RCI (7 days w/o RCI)	Daily	10494040_EN Rev. B, 2011-02
TDM	Tobramycin [TOB_2]	10376932 [10376932]	400 [4 x 100]	S, P (EDTA, Hep Citrate, Ox/Fl)	0.4–13.5 µg/mL [0.9–28.9 µmol/L]	Enzyme multiplied immunoassay technique (EMIT)	30 days	Drug Cal 2 10376771	30 days	Daily	10494043_EN Rev. B, 2011-02
TDM	Valproic Acid [VPA_2]	10377508 [10377508]	400 [4 x 100]	S, P (EDTA, Hep Citrate, Ox/Fl)	3.0–176.4 µg/mL [20.8–1222.5 µmol/L]	Enzyme multiplied immunoassay technique (EMIT)	30 days with RCI [14 days w/o RCI]	Drug Cal 2 10376771	30 days with RCI (14 days w/o RCI)	Daily	10494059_EN Rev. B, 2011-02
TDM	Vancomycin [VANC_2]	10377509 [10377509]	400 [4 x 100]	S, P (EDTA, Hep Citrate, Ox/Fl)	1.0–52.6 µg/mL [0.7–35.2) µmol/L]	Enzyme multiplied immunoassay technique (EMIT)	30 days	Drug Cal 2 10376771	30 days	Daily	10494057_EN Rev. B, 2011-02
Toxicology	Acetaminophen [ACET]	07989138 [10327381]	600 [4 x 150]	S, P (Li Hep)	0.2–20 mg/dL [13.2–1322 µmol/L]	Enzymatic	28 days with RCI [10 days w/o RCI]	Tox/Ammonia Cal 10309217	28 days with RCI (10 days w/o RCI)	Daily	10493934_EN Rev. F, 2014-06
Toxicology	Ethanol [ETOH_2]	10379635 [10379635]	600 [4 x 150]	S, P (EDTA, Li Hep)	10.0–600.0 mg/dL [2.2–130.2 mmol/L]	Enzymatic	30 days	Tox/Ammonia Cal 10309217	14 days	Daily	10493987_EN Rev. B, 2011-06
Toxicology	Salicylate [SAL]	07989456 [10327382]	600 [4 x 150]	S, P (Li Hep)	3.0–100.0 mg/dL [3.0–100.0 mg/dL]	Enzymatic	60 days	Tox/Ammonia Cal 10309217	60 days	At calibration	10494029_EN Rev. E, 2011-02
Toxicology	Serum Barbiturates [sBRB]	05040084 [10335934]	400 [4 x 100]	S, P (Li Hep)	0.06–6.0 µg/mL [0.06–6.0 µg/mL]	Homogenous competitive enzyme immunoassay	30 days	Toxicology Cal 1 10283016 Toxicology Cal 2 10283017 Toxicology Cal 3 10283018 Toxicology Cal 4 10283019	7 days	At calibration	10494031_EN Rev. B, 2011-02
Toxicology	Serum Benzodiazepines [sBNZ]	05042192 [10335939]	400 [4 x 100]	S, P (Li Hep)	10–200 ng/mL [10–200 ng/mL]	Homogenous competitive enzyme immunoassay	30 days	Toxicology Cal 1 10283016 Toxicology Cal 2 10283017 Toxicology Cal 3 10283018 Toxicology Cal 4 10283019	7 days	At calibration	10494030_EN Rev. B, 2011-02
Toxicology	Tricyclic Anti-depressants [TCA]	06520218 [10283014]	400 [4 x 100]	S, P U	80–1000 ng/mL [80–1000 ng/mL]	Homogenous competitive enzyme immunoassay	30 days	Toxicology Cal 1 10283016 Toxicology Cal 2 10283017 Toxicology Cal 3 10283018 Toxicology Cal 4 10283019	7 days	At calibration	10494038_EN Rev. B, 2011–01

P: Plasma
RCI: Reagent Container Insert

S: Serum
U: Urine

ADVIA, CardioPhase, Emit, and all associated marks are trademarks of Siemens Healthcare Diagnostics Inc. or its affiliates. All other trademarks and brands are the property of their respective owners.

Product availability may vary from country to country and is subject to varying regulatory requirements. Please contact your local representative for availability.

Order No. A91DX-CAI-150023-GC1-4A00
04-2015 | All rights reserved
© 2015 Siemens Healthcare Diagnostics Inc.

Global Siemens Headquarters

Siemens AG
Wittelsbacherplatz 2
80333 Muenchen
Germany

**Global Siemens Healthcare
Headquarters**

Siemens AG, Healthcare
Henkestrasse 127
91052 Erlangen
Germany
Telephone: +49 9131 84-0
www.siemens.com/healthcare

Global Division

Siemens Healthcare Diagnostics Inc.
511 Benedict Avenue
Tarrytown, NY 10591-5005
USA
www.siemens.com/diagnostics